

Electrical Lice Remover

Presented by KaniStar Group

INTRODUCTION

Electrical Lice Remover is an advanced device for the treatment of head lice and nits. It is developed by KaniStar Group and allows you to combine the power of blowing with a lice comb and the application of controlled heat for the complete elimination of lice without damaging the hair, all in one machine!

This combination allows your business to be more efficient by needing **only one machine for the whole treatment**, with the added benefit of being a totally natural, non-chemical way of eliminating head lice and nits.

The main heating unit offers a multitude of beneficial functions and features.

1. The operating noise level is very low. (53dB).
2. It takes less than 3.5 min for the temperature to reach 60 °C .
3. Easy to handle and maintain, smart trouble shooting system installed.
4. Accurate LCD display, real time temperature of comb monitoring.
5. The unit can be freestanding.
6. Compact and light weighted, using little space, easy to transport.
7. Swiveling easy to fit outlet port, makes attaching the hose effortless (Both the removal and attachment should be done by professional persons.).
8. The timer will sound an alarm each 35 seconds to remind the user for preventing overheating temperature on the scalp.

Package Inspection

The day the unit is delivered, an inspection of the outer packing should be made first for any damage to the packing box.

Then the unit should be unpacked carefully, all components and the machine should be inspected for any damage.

If any damage is evident, the carrier or dealer should be notified immediately.


The manufacturer is not liable for any damage caused by a third party "transport or delivery company".

Operation Instructions

The unit should be placed in a well ventilated position and all inlets and outlets should be unobstructed. There should be a 10cm or 4-inch clearance on all sides.

Connect the air flow hose to the air outlet.

Connect the power plug to a suitable power outlet, turn on the switch on the back and the unit is now in standby mode.

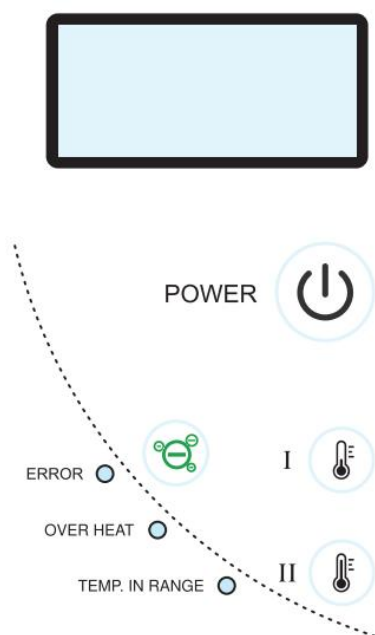

Press and hold the Power button on the control panel. The unit will begin to blow heated air; and countdown will be displayed on the LCD screen.

Press the button I/II will adjust the heating power in 2 levels. Level I is 55°C. Level II is 59°C.

Heating can be cancelled by pressing the same start button and the unit will be shut down in 2 mins.

The timer will sound an alarm each 35 seconds to remind the user for preventing overheating temperature on the scalp.

Its heat is controlled to eliminate lice and nits.

SAFETY WARNING

All maintenance and repairs should be carried out by qualified authorized trained technicians.

All parts' replacement should be done or supervised by a qualified trained technicians

If the unit starts to make a continuous noise, please contact an Authorized service centre.

The machine housing can be wiped clean with a mild cleaning agent that is non-abrasive.

SPECIFICATIONS

Dimension 30cm×22cm×36cm (12" ×8.65" ×14")

Noise Level 49dB (high speed)

43dB (low speed)

Negative ion $6 \times 10^6 \text{Pcs/cm}^3 \pm 10\%$

Weight 6.14kg (13.52lbs)

Vent Hose Food class, flexible, sensor inside, 1.5m long.

Suggested filter replacement cycle 2 months or over 120 hours

Suggested Operation Environment

Temperature: -20°C -30°C

Relative Humidity 35% -80%

Altitude $\leq 3200\text{m}$

Warming temperature: 55/59°C

Warming time Less than 4 minutes

Alarming temperature: 65°C

Temperature Control Tolerance $\pm 1.0^\circ\text{C}$

Storage and Transportation Temperature -10°C -45°C

ELECTRICAL PARAMETERS

Power Source: 110V, 60Hz/220V, 50Hz

Heating Element PTC, 650W (110V) /700W (220V)

Fan Motor	45W
Fuse	Porcelain Fuse
Regulation	IEC60601
Rotational Speed	110V / max 18m/s 220V / max 18m/s
Rated Fuse	15A (125VAC) 8A (250VAC)

TROUBLE SHOOTING

Overheat

Once temperature exceeds $\geq 65^{\circ}\text{C}$, the red alarm light will flash and the alarm will sound, heating elements and fans will be shut down automatically, the operation panel is not workable.

Unplug the power plug to allow sufficient cooling time before attempting to restart.

If the problem persists, please contact your KaniStar Authorized service centers.

Fault Alarm

The yellow light flashing and alarm sounding indicate possible malfunctions with mechanical system, fans, heater or sensor, etc.

Temperature should be closely monitored until the machine can be turned off and unplugged.

Wait ten minutes and restart, the machine will run self-diagnosis; if the fault light remains on, please consult your local distributor or service centre.

Power Failure

Check that the power cable is connected.

Check that the Power supply is connected.

Check the fuse.

PACKING LIST

Air-conditioning Unit

Instruction Manual

Rubber Tubing Set

Power Cord

WARRANTY

12 months limited warranty from the date of purchase.

Warranty does not cover any misuse of the item or any damage caused directly or indirectly by the user.

An extended warranty can be purchased from the distributor.

Dallas

LLC. 3249 Story Rd. W Irving, TX 75038

Tel: +1 469 930 6359

Web: www.aeoluspet.com

Shanghai

Floor 5, Building 6, No. 675, Zhujin Road, Songjiang

District, Shanghai, China.

Tel: +86-21-67632060, 67632364

Web: www.toexpet.com